

TN_P10; Anna Gakamsky

Introduction

The absorbance A or optical density (OD) of a substance is inherent in all aspects of spectroscopy and follows the Beer-Lambert law¹, $A(\lambda) = \log[T_{ref}(\lambda)/T_s(\lambda)]$, where $T_{ref}(\lambda)$ is the transmission of the reference sample and $T_s(\lambda)$ the transmission of the sample at the measured wavelength.

Instrument and procedure

The FS5 Fluorescence Spectrometer is equipped with a silicon detector that allows for transmission and absorbance measurements. Liquid samples can be measured in the standard cuvette holder (SC-05), cuvette holder for coolant circulation (SC-20) and thermoelectrically cooled/heated cuvette holder (SC-25) which can be interchangeably installed in the FS5. Fluorescein solutions were prepared in phosphate buffered saline (PBS) in 1 cm plastic cuvettes. After the transmission scans of the sample and reference are complete, they can be joined to calculate the absorbance through the absorption wizard in the operating software, Fluoracle.

Results

Figure 1 shows the transmission spectra of fluorescein samples at different concentrations measured in an FS5 fluorescence spectrometer.

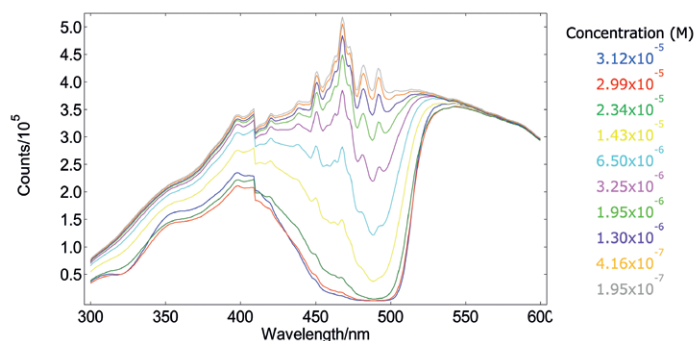


Figure 1: Transmission spectra for different concentrations of fluorescein in PBS.

The transmission spectra were recorded with a band-pass of 3 nm, a step of 1 nm and integration time of 0.2 s. The repeatability of the measurements was confirmed by reference transmission scans to less than 1% relative deviation.

From the transmission spectra, the absorbance spectra of the fluorescein samples in PBS shown in Figure 2 were calculated in the Fluoracle software. The spectra reveal the dianion form² of the solution with the absorption peak at 490 nm. The concentration c was calculated according to the Beer-Lambert law from the peak absorbance as $c = A / \epsilon l$, with the molar attenuation coefficient ϵ of fluorescein².

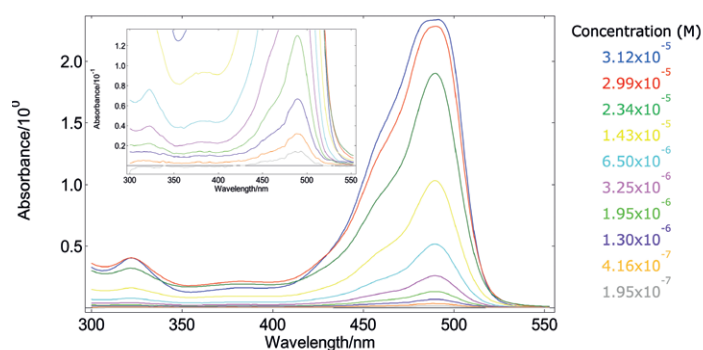


Figure 2: Absorbance spectra for different concentrations of fluorescein in PBS.

Conclusion

The FS5 Fluorescence Spectrometer enables transmission and absorbance measurements of liquid samples in the wavelength range 230 nm-1000 nm and optical densities of OD 0.01 to OD 2.

TN_P10; Anna Gakamsky

References

1. Ball, D. W. *The Basics of Spectroscopy*. (SPIE Press, 2001).
2. Sjöback, R., Nygren, J. & Kubista, M. Absorption and fluorescence properties of fluorescein. *Spectrochimica Acta Part A: Molecular and Biomolecular Spectroscopy* **51**, L7–L21 (1995).

Proprietary notice

Words and logos marked with © or ™ are registered trademarks or trademarks owned by Edinburgh Instruments Limited. Other brands and names mentioned herein may be the trademarks of their respective owners.

Neither the whole nor any part of the information contained in, or the product described in, this document may be adapted or reproduced in any material form except with the prior written permission of the copyright holder.

The product described in this document is subject to continuous developments and improvements. All particulars of the product and its use contained in this document are given by Edinburgh Instruments in good faith. However, all warranties implied or expressed, including but not limited to implied warranties of merchantability, or fitness for purpose, are excluded.

This document is intended only to assist the reader in the use of the product. Edinburgh Instruments Limited shall not be liable for any loss or damage arising from the use of any information in this document, or any error or omission in such information, or any incorrect use of the product.

Confidentiality status

This document is Open Access. This document has no restriction on distribution.

Feedback on this Application Note

If you have any comments on this Application Note, please send email to sales@edinst.com giving:

the document title

the document number

the page number(s) to which your comments refer

an explanation of your comments.

General suggestions for additions and improvements are also welcome.

For more information, please contact:

T: +44 (0) 1506 425 300 E: sales@edinst.com

F: +44 (0) 1506 425 320 W: www.edinst.com